-online.com genomics



2 Images

|--|

Quantity:	500 units
Application:	Polymerase Chain Reaction (PCR)

Product Details

Purpose:	TransFast® Taq DNA Polymerase is an engineered version of Taq DNA Polymerase.
Brand:	TransFast®
Specificity:	The enzyme consists of a single polypeptide with a molecular weight of approximately 94 kDa. TransFast® Taq DNA Polymerase has 5'-3' DNA polymerase activity and 5'-3' exonuclease activity. It lacks 3'-5' exonuclease activity.
Characteristics:	- Extension rate is about 6 kb/min. - Template-independent "A" can be generated at the 3' end of the PCR product. PCR products can be directly cloned into pEASY®-T vectors. - Amplification of genomic DNA fragment up to 4 kb.
Components:	DNA Polymerase, 10X Taq Buffer, 6X DNA Loading Buffer
Unit Definition:	One unit of TransFast® Taq DNA Polymerase incorporates 10 nmol of deoxyribonucleotide into acid-precipitable material in 30 minutes at 74°C.

Application Details

Application Notes:	Routine PCR, High throughput PCR,Colony PCR
Comment:	TransFast® Taq DNA Polymerase has passed the following quality control assays: functional
	absence of double- and single-strand endonuclease activity, >99% homogeneous measured by
	SDS-PAGE. Each batch of TransFast $\ensuremath{^{ m B}}$ Taq DNA Polymerase has been assayed for
	amplification efficiency to amplify p53 gene from 10 ng of human genomic DNA.

Order at www.genomics-online.com USA & Canada: +1 877 302 8632 | support@antibodies-online.com Page 1/3 | Product datasheet for ABIN5519490 | 09/12/2023 | Copyright antibodies-online. All rights reserved.

Application Details	
Restrictions:	For Research Use only
Handling	
Buffer:	Storage Buffer: 20 mM Tris-HCl (pH 8.0), 0.1 mM EDTA, 1 mM DTT, 100 mM KCl, 50 % glycerol,
	stabilizers.
	10xTransFast® Taq Buffer (with Mg2+): 200 mM Tris-HCl (pH 8.4), 100 mM KCl, 100 mM
	(NH4)2SO4, 20 mM MgSO4, others
Storage:	-20 °C
Storage Comment:	at -20°C for two years
Expiry Date:	24 months
Publications	
Product cited in:	Johnson, Drugan, Miller, Evans: "38" in: , Vol. 1363, Issue Nucleic acids research, pp. 28-39, (
	1991)

Images



M: 1Kb Plus DNA Ladder E: *EasyTaq®* DNA Polymerase F: *TransFast® Taq* DNA Polymerase Image 1.

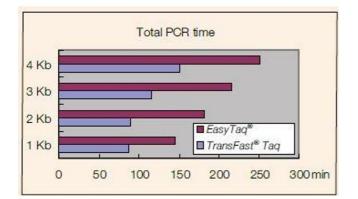


Image 2.